

↑ TO REMOVE - PULL HERE ↑

Utah Newborn Screening


Please Read This Important Information

MISCELLANEOUS:

First screen specimens should be collected on ALL newborns between 48 hours and 5 days of age. If newborn is discharged before 48 hours obtain specimen within 4 hours of discharge. A second screen is required between 7 days and 28 days of age.

COLLECTION INSTRUCTIONS

1. Legibly print ALL information in spaces provided using block capital letters.
2. Mark "test requested" box to indicate which screen (first or second) needs to be run.
3. If using this form to replace a lost or unusable kit form, cross out the miscellaneous number and write in original kit number.
4. Collect specimen with heel stick. See newborn screening handbook for detailed instructions.
5. Fill all 7 circles.
6. Dry 3-4 hours before mailing.

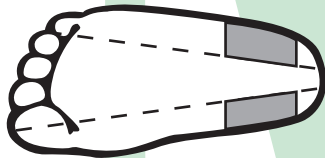

COLLECT SAMPLE FROM
SHADED AREA.

Mailing Instructions:

1. When blood is dry, fold stock card (from back of form) over blood spots. The flap should enclose the blood spots and reveal a biohazard symbol.
2. If using the postal service, place form with blood spots covered into envelope.

Mail to: Newborn Screening Laboratory
Utah Department of Health
46 N Medical Dr
Salt Lake City UT 84113-9903
Phone: (801) 584-8256

For more information, call, refer to your
handbook or visit our website:

http://www.health.utah.gov/newbornscreening/HCP_Instructions.htm

Retain this sheet for your records.

I.D. Number